

**Summary of example:**

In this example we estimate some basic synaptic parameters based on a simulated voltage clamp experiment with synaptic input.

The target data (iclamp.dat) consist of the recorded clamp current (in units of nA) from a virtual voltage clamp electrode inserted into a single-compartment model, which contains Hodgkin-Huxley- type Na<sup>+</sup>, K<sup>+</sup> and leak conductances, plus a conductance-based synapse with a double exponential time course (rise time: 0.3 ms, decay time: 3 ms, maximal conductance: 10 nS, delay: 0 ms). The data file also contains time.

The model neuron receives through the synapse a spike train input, which consists of 4 spikes at regular 100 ms intervals starting at 100 ms. The full length of the recording is 500 ms, the sampling rate is 40 kHz.

The model file (simple\_hh\_syn\_vclamp\_toopt2.hoc) contains the neuronal model (including the synapse), the spike generator (NetStim) object which generates the input, and a NetCon object which connects the input to the cell.

As we need to set the parameters of the synapse and those of the NetCon, and Optimizer cannot discover these parameters automatically, we use a simple user function (ufun.txt) to adjust the parameters (maximal conductance (in microsiemens), synaptic delay, rise time constant, decay time constant (all in milliseconds)).

We need voltage clamp at a constant level (-70 mV); one way to accomplish this is to use a step protocol in voltage clamp with a single amplitude of -70 mV (and arbitrary delay and duration), and an initial voltage of -70 mV.

The optimization should be done using the mean square error cost function.

**input file:** iclamp\_new\_evendtdat

**model:** simple\_hh\_syn\_vclamp\_toopt2.hoc

**Step-by-step instructions to run the example from the Neuroptimus GUI:**

Run „python3 neuroptimus.py -g” to start the GUI

Neuroptimus

Menu

Target data

Model

Settings

Fitness

Run

Results

Statistics

Input Type

Current trace

Data file

irrogate/iclamp\_new\_eventdt.dat

Browse...

☒ Contains time

Base directory

neuroptimus/new\_test\_files/testcase\_2\_VClamp\_surrogate

Browse...

Number of traces

1

Units

nA

Length of traces (ms)

500


Sampling frequency (Hz)

40000.0

Load data


Current traceiclamp\_new\_eventdt.dat

At 'Data File' load the target data, at 'Base Directory' choose the directory where you want to save the results. Fill out all the cells and press 'Load data'.Go on by pressing the Model tab.

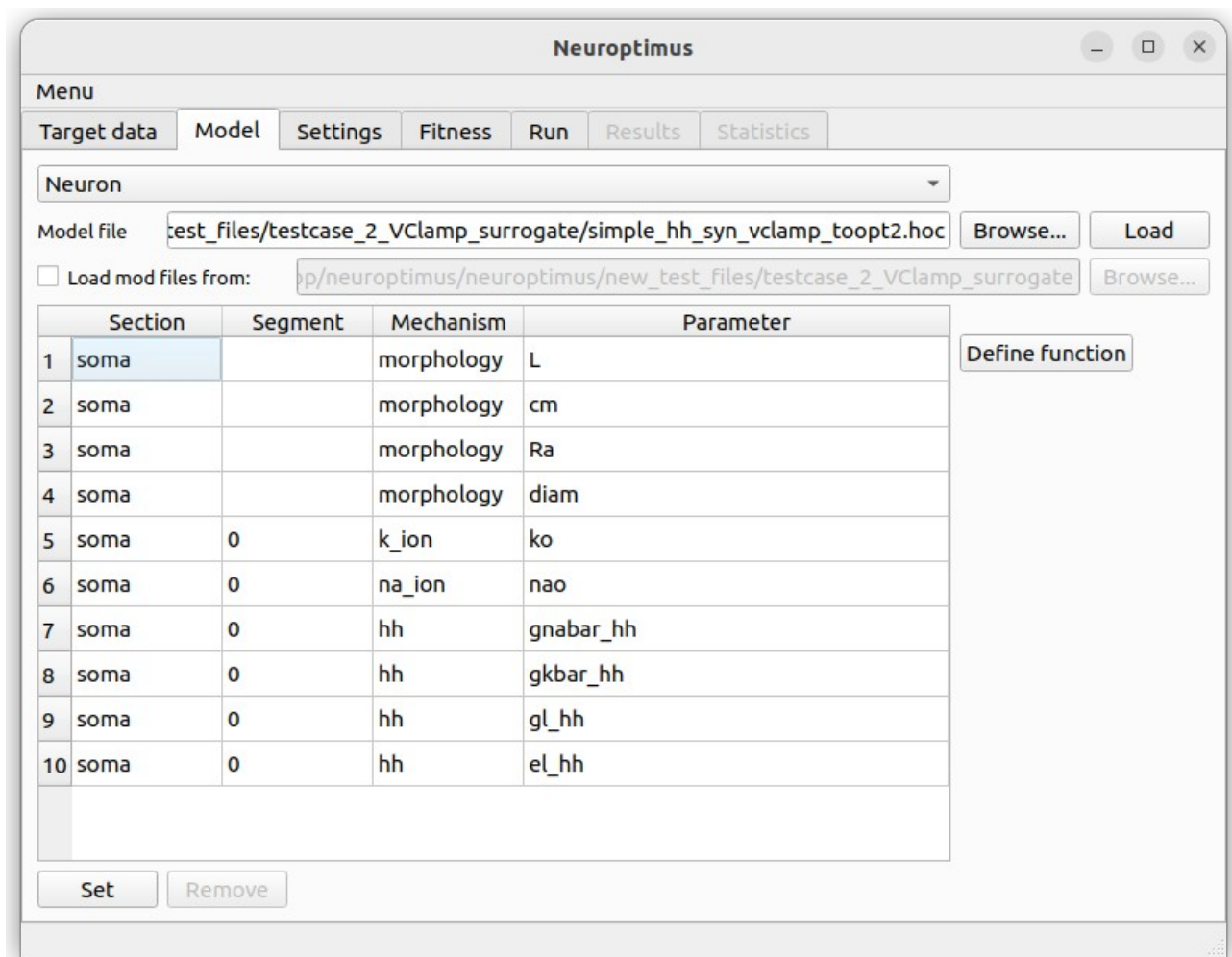

Browse to the model file and load the model.

Press the "Define function" button. In the new window, press 'Load' to select the user-defined function 'ufun.txt'.

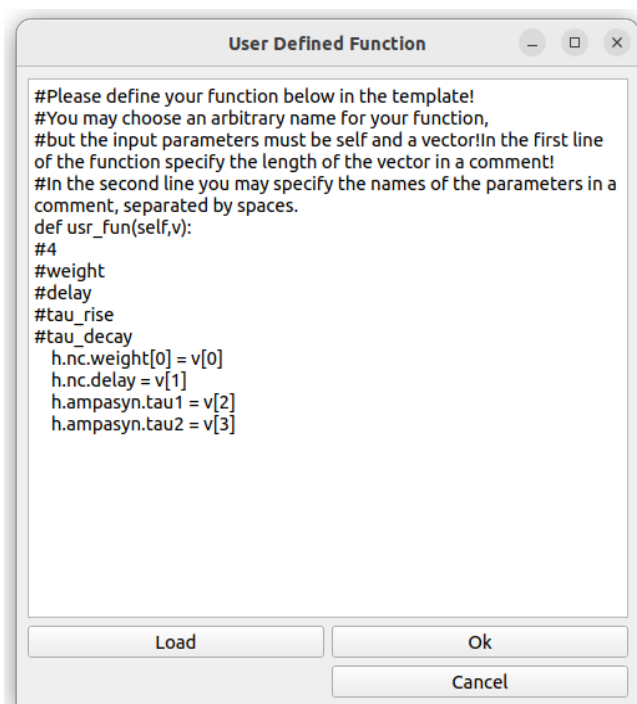

Press 'Ok', then go on by pressing the Settings tab.

Neuroptimus

Menu

Target data

Model

Settings

Fitness

Run

Results

Statistics

Stimulation protocol

VClamp

Stimulus Type

Step Protocol

Amplitude(s)

Stimulation Time Settings

Delay (ms)

1

Duration (ms)

1

Stimulus Position Configuration

Section

soma

Position inside section

0.5

Recording Settings

Parameter to record

i

Section

soma

Position inside section

0.5

Simulation Settings

Initial voltage (mV)

-70

tstop (ms)

500.0

Time step

0.025

Fill in all the cells. Press 'Amplitude(s)' to open a new window. and set the amplitude of the stimulus.

Stimuli Window

Number of stimuli:

1


Create

Amplitude (mV)


1	-70
---	-----

Accept

Go on by pressing the Fitness tab.

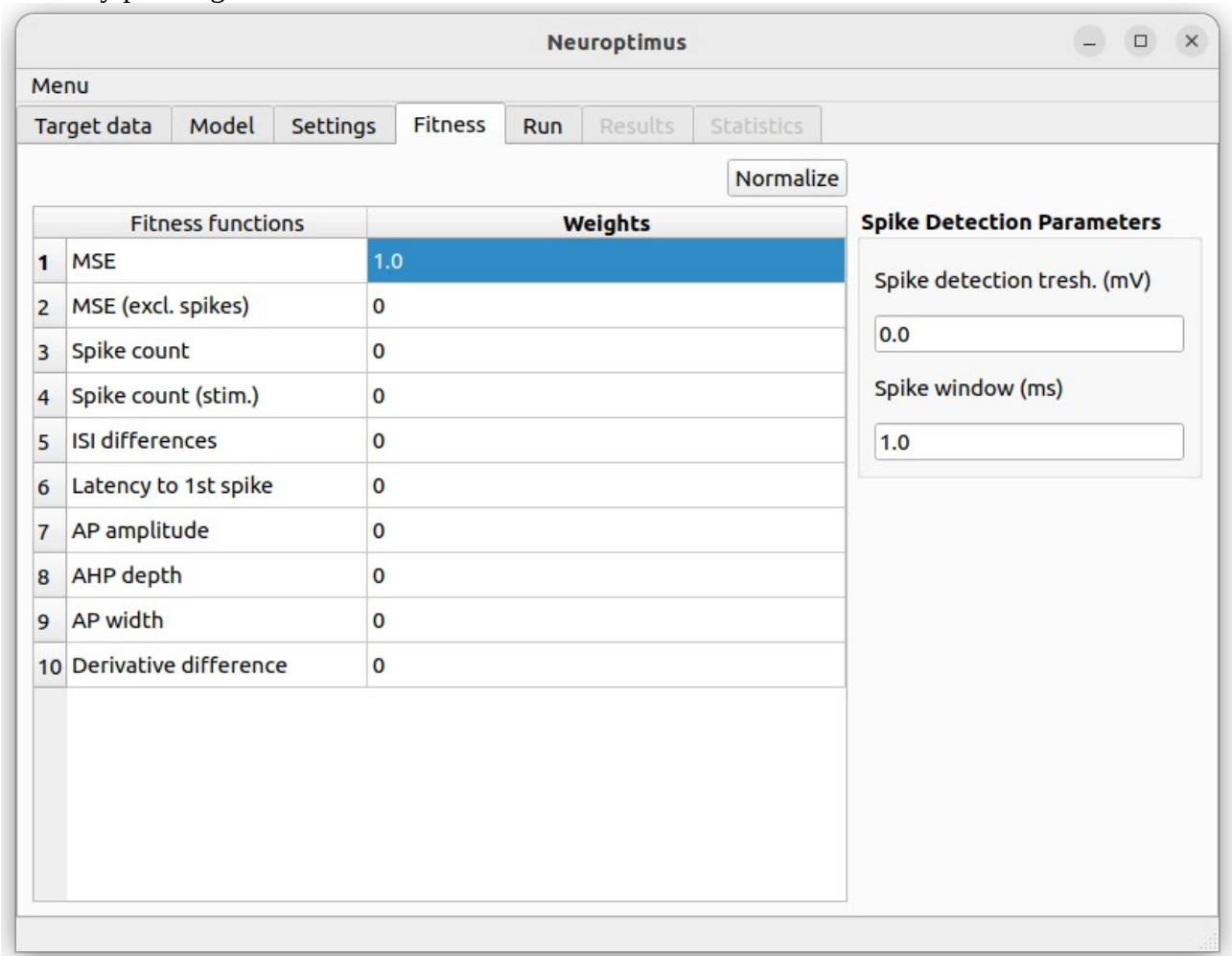

Choose the 'MSE' fitness function, and set its weight to 1. Go on by pressing the Run tab.

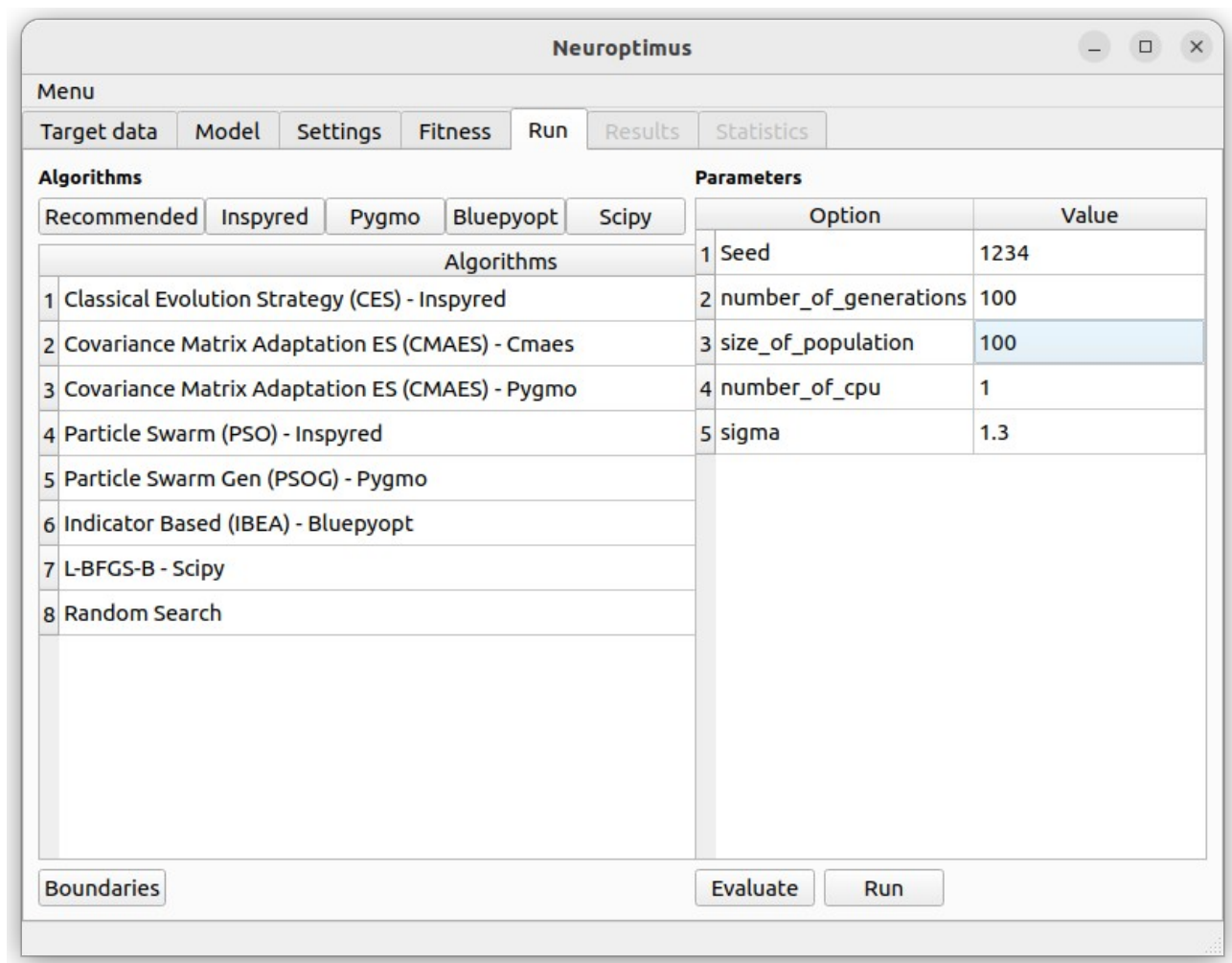


Select an algorithm, and set its parameters.

Press the 'Boundaries' button to define the boundaries of the parameters to be optimized:

Press 'Set'.

Boundaries can also be loaded from a file.

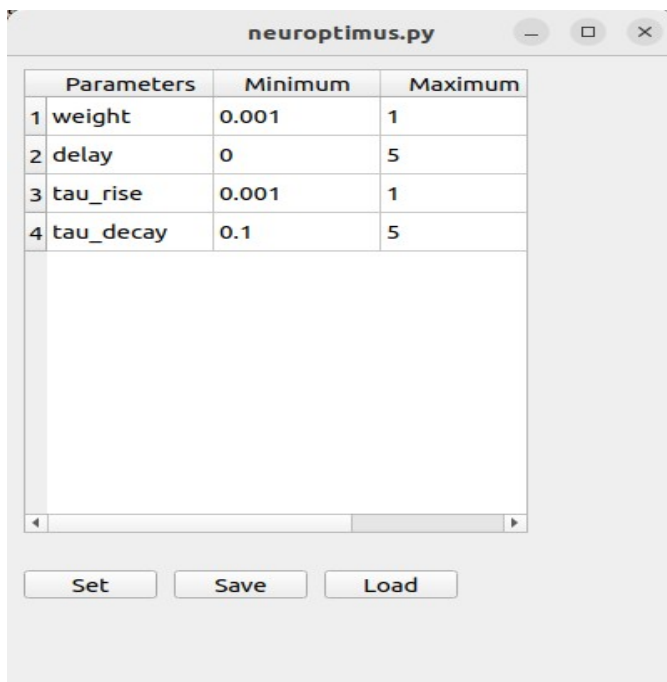

Start the optimization pressing the 'Run' button.